


Capsule Hooks

Capsule hooks are used for temporary support of the capsular bag during a cataract surgery when the zonules are weak or in case of subluxated lenses.

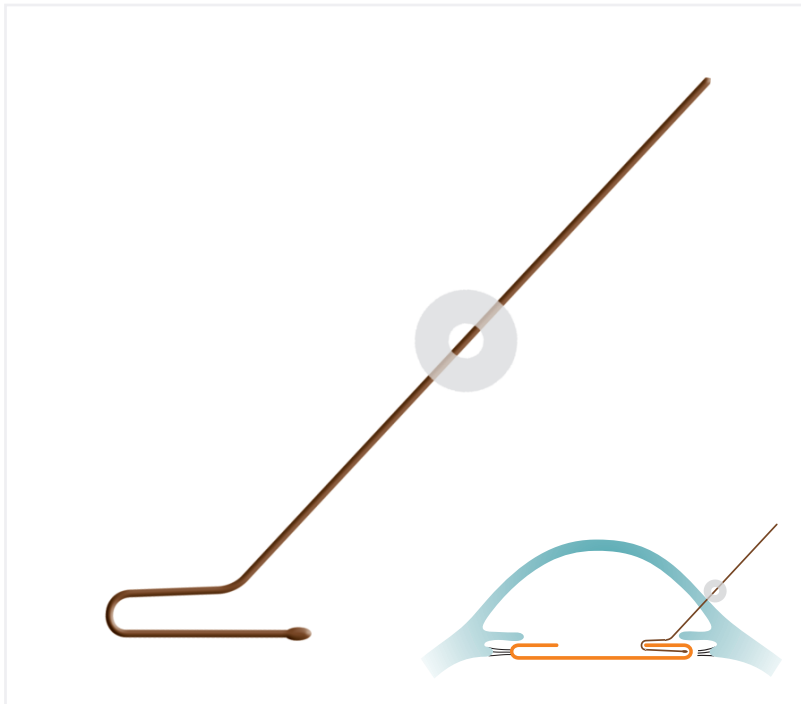

Product Description

Product Code	Description	Qty. per box
MIPL/D8	Capsule Hooks	5 Pcs. (Sterile)

Features

- Capsule hooks are made of Nylon with an adjustable silicon lock
- Designed such that the hook is angled to be in the plane of the anterior capsule
- Easy to get a firm and reliable grip on the rim
- Tip of the hook is rounded to prevent any accidental damage to the peripheral capsule while placing the hooks in the bag
- The round tip also provides a larger support area per hook there by reducing the number of hooks required in each surgery
- Silicon lock helps in retaining the hook in the retracted position during the surgery
- CE certified

Advantages

- Maintains the position of the capsule during phaco
- Supports the fornix as well as the capsulorhexis rim, so that the peripheral capsule isn't attracted towards the phaco tip
- Avoids the problem of the capsular tension ring trapping cortex
- Prevents zonular stress during nuclear rotation
- Makes it possible to perform posterior chamber phaco even with subluxated lenses having minimum zonular support
- Easy to insert and remove

Instructions for Use

The instructions for use and surgical technique are supplied in the form of a leaflet with instructions given in diagram form where appropriate . The international symbols used on the packaging are also explained.

Sterilization

Products are sterilized through Ethylene Oxide (EO) according to a validated cycle.



MADHU INSTRUMENTS PVT LTD

A-260, Okhla Industrial Area, Phase - 1, New Delhi - 110020, India

T: (+91-11) 4270 1030, 31 | E: sales@madhuinstruments.com | www.madhuinstruments.com

Customer Care No. +91 9810866176 (for international queries), +91 7428186111 (for domestic queries)