


Vitreotomy Infusion Cannula

Disposable cannula used for infusion of irrigation fluids or surgical fluids like silicon oil/PFCL during vitreo-retina ophthalmic surgical procedures.

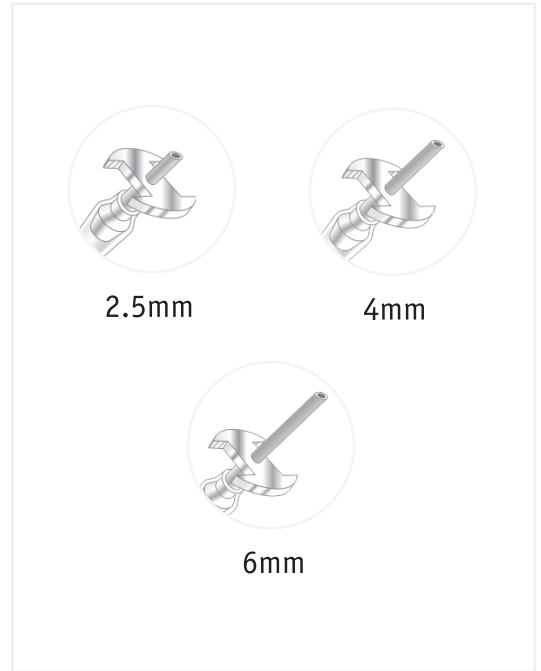

Product Description

Product Code	Description	Qty. per box
MIPL/E1/2.5	2.5 mm tip (20G)	10 pcs. (Sterile)
MIPL/E1/4.0	4.0 mm tip (20G)	10 pcs. (Sterile)
MIPL/E1/6.0	6.0 mm tip (20G)	10 pcs. (Sterile)

Features

- Thin wall stainless steel tube is used for maximum flow
- Designed especially for silicon oil injection
- Superior grip of silicon tube at cannula end as well as connector end
- CE certified

Instructions for Use

The instructions for use and surgical technique are supplied in the form of a leaflet with instructions given in diagram form where appropriate . The international symbols used on the packaging are also explained.

Sterilization

Products are sterilized through Ethylene Oxide (EO) according to a validated cycle.



MADHU INSTRUMENTS PVT LTD

A-260, Okhla Industrial Area, Phase - 1, New Delhi - 110020, India

T: (+91-11) 4270 1030, 31 | E: sales@madhuinstruments.com | www.madhuinstruments.com

Customer Care No. +91 9810866176 (for international queries), +91 7428186111 (for domestic queries)